



DAS 8016

Pocket Reader System

USB Charging &
USB Serial Data
Interface



Bluetooth Wireless
Keyboard Output

Introducing our new DAS-8016 pocket reader system with wireless capability. The reader can scan and display transponder ID on the bright OLED display. At 1/2 the size of our current reader systems, the DAS-8016 provides a fast and pocketable reader for today's mobile researcher.

All data scanned can be stored or transmitted directly to nearby devices using several output options; Direct USB serial with cable for secure communication, Bluetooth HID keyboard connection for ease of connection to any device, or paired with a BMDS Wireless COM Module for maximum flexibility of both Keyboard and Serial USB connections. The DAS-8016 is compatible with our IMI family of Read-Only transponders (IMI-400, IMI-500, and IMI-1000)

Features:

- Compatible with IMI-400 (6mm), IMI-500 (8mm), and IMI-100 (11mm)
- Serial USB port output
- Bluetooth Wireless HID Keyboard output
- Compatible with BMDS Wireless Com Module
- Memory storage for 500 records
- Audio feedback
- USB charging for maximum flexibility




DAS

8016

Pocket Reader System

Optional Accessories:

- **BMDS Wireless Com Module** – Optional wireless communication with both USB Serial and USB Keyboard Output

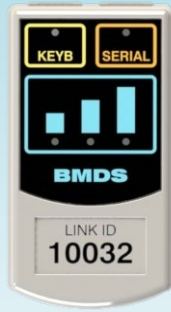

Reader Specifications:

- Read Only – Display of ID from IMI-400, IMI-500, and IMI-1000 transponders
- Super Bright OLED display
- Data Storage and Download of 500 IDs
- Onboard Menu Driven Setup
- USB Port for secure serial connection and recharging of Lithium Ion battery
- Selectable Record Format Parameters

- Bluetooth transmission distance approximately 20' (approximately 6 m)
- Wireless Com Module transmission over secure P2P Protocol 903MHz ISM Band. Receive distance approximately 20' (approximately 6 m)
- Transponder read distance 2" (approximately 6 cm)
- 1.0 A USB Power Supply for Recharging
- Battery: Single 3.7 Li-ion 3000mAH User Replaceable
- Reader: 2.1" Wide by 1.6" High by 8.5" Long

Reference

To the right, is a size comparison of the current full feature portable Wireless Reader and the new simpler Pocket Reader.



CALL US: 800 526 BMDS (2637)