




DAS 8017

Pocket Reader System

USB Charging &
USB Serial Data
Interface



Bluetooth Wireless
Keyboard Output

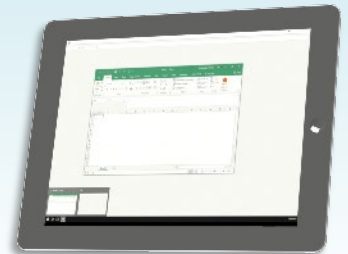

Introducing our new DAS-8017 pocket reader system with wireless capability. The reader can scan and display transponder ID on the bright OLED display. At 1/2 the size of our current reader systems, the DAS-8017 provides a fast and pocketable reader for today's mobile researcher.

All data scanned can be stored or transmitted directly to nearby devices using several output options; Direct USB serial with cable for secure communication, Bluetooth HID keyboard connection for ease of connection to any device, or paired with a BMDS Wireless COM Module for maximum flexibility of both Keyboard and Serial USB connections. The DAS-8017 is compatible with our IPT-300, IPTT-300, and XPT-100 family of programmable transponders.

Features:

- Compatible with XPT-100 (8mm) and IPT/IPTT-300 (14mm) transponders. Does not display IPTT Temperature Data
- Serial USB port output
- Bluetooth Wireless HID Keyboard output
- Compatible with BMDS Wireless Com Module
- Memory storage for 500 records
- Audio feedback
- USB charging for maximum flexibility




DAS

8017

Pocket Reader System

Optional Accessories:

- **BMDS Wireless Com Module** – Optional wireless communication with both USB Serial and USB Keyboard Output

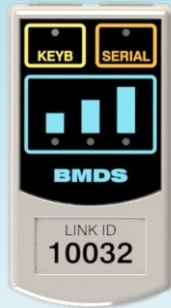

Reader Specifications:

- Read Only – Display of ID from XPT-100, IPT-300, and IPTT-300 transponders
- Super Bright OLED display
- Data Storage and Download of 500 IDs
- Onboard Menu Driven Setup
- USB Port for secure serial connection and recharging of Lithium Ion battery
- Selectable Record Format Parameters

- Bluetooth transmission distance approximately 20' (approximately 6 m)
- Wireless Com Module transmission over secure P2P Protocol 903MHz ISM Band. Receive distance approximately 20' (approximately 6 m)
- Transponder read distance 2" (approximately 6 cm)
- 1.0 A USB Power Supply for Recharging
- Battery: Single 3.7 Li-ion 3000mAH User Replaceable
- Reader: 2.1" Wide by 1.6" High by 8.5" Long

Reference

To the right, is a size comparison of the current full feature portable Wireless Reader and the new simpler Pocket Reader.



CALL US: 800 526 BMDS (2637)